

OTI

**OTIScan 3000
Installation Instructions**


Version 1.0
12 September 2007

Contents


OTIScan 3000 Application Installation	2
OTIScan 3000 Driver Installation	4
AScan Driver Installation	8
Copy msvc71.dll from CD	12
OTIScan 3000 Uninstallation	13

OTIScan 3000 Application Installation


1. Insert the OTIScan 3000 CD and double click on the OTIScan 3000 installer to run it.
2. Click the Next button on the Welcome screen.

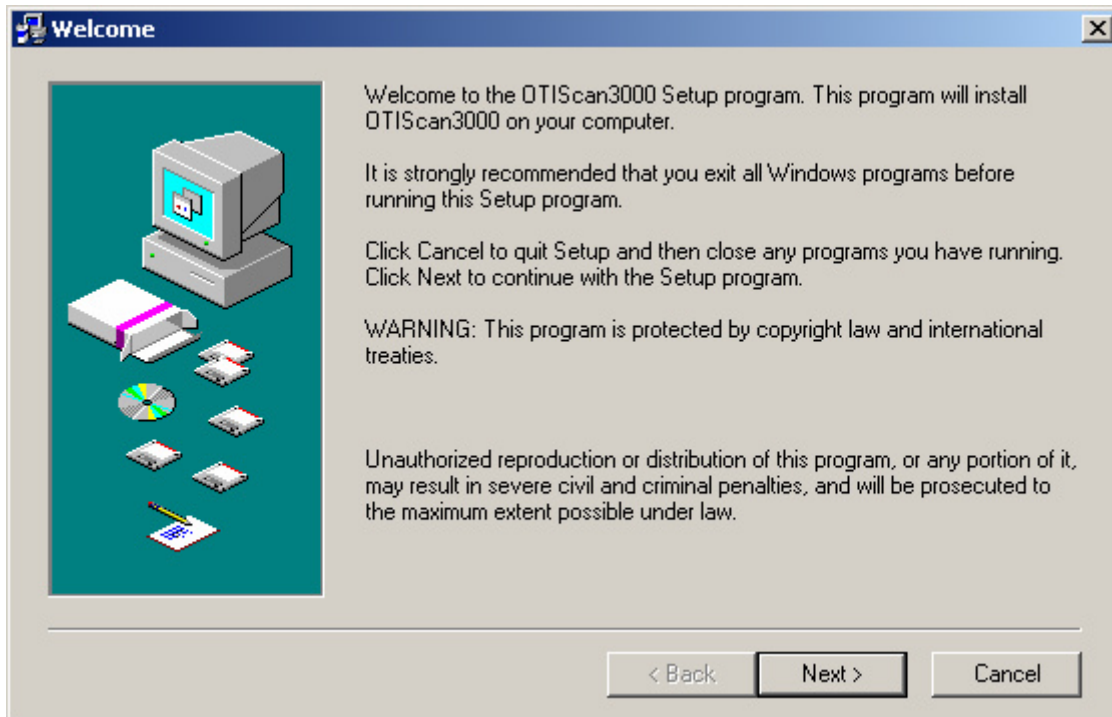


3. If templates from a previous OTIScan 3000 installation are detected, the Template Options screen will be displayed. Leave the checkbox on this screen checked unless you have modified your templates and do not wish to overwrite them. Click the Next button.

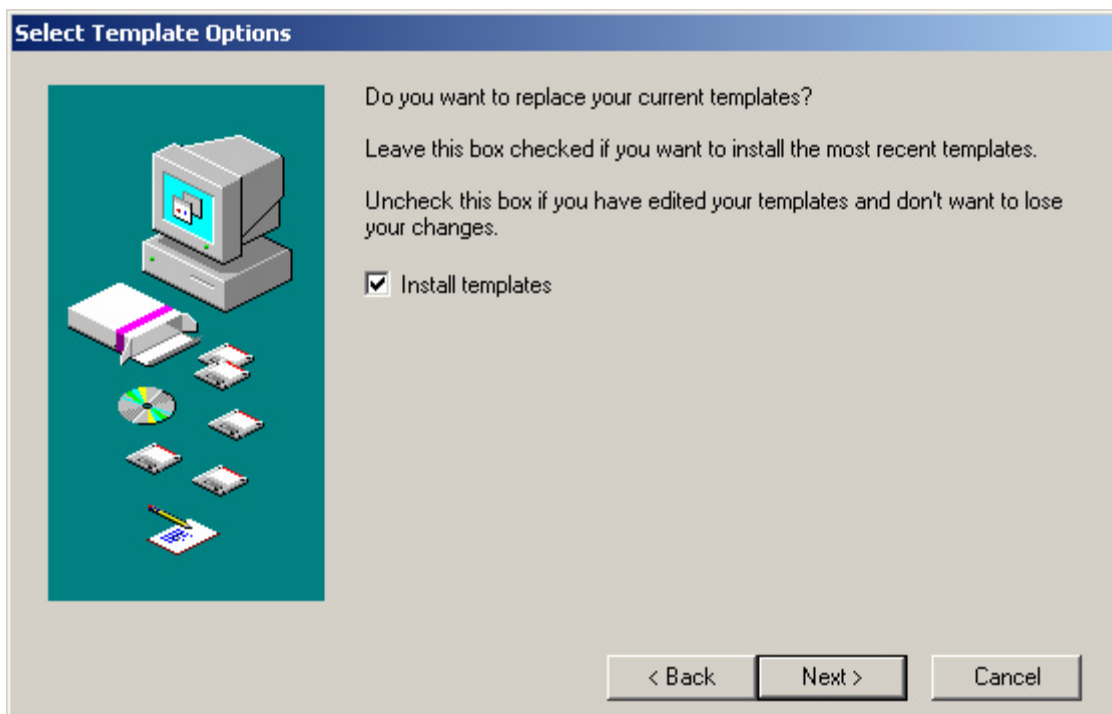

4. Click on the Next button on the Ready To Install screen. OTIScan 3000 will be installed.

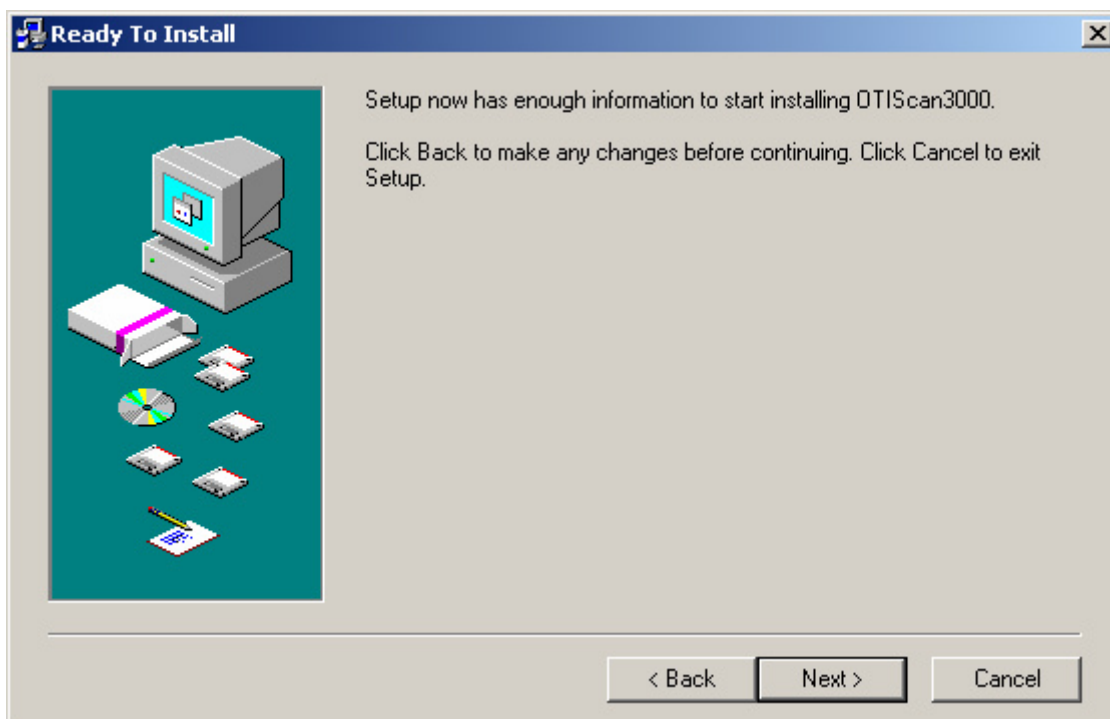


5. On the Finished screen is a checkbox to install the OTIScan 3000 device driver. Uncheck it if you have previously installed the OTIScan 3000 driver and do not wish to re-install it. Click on the Close button.

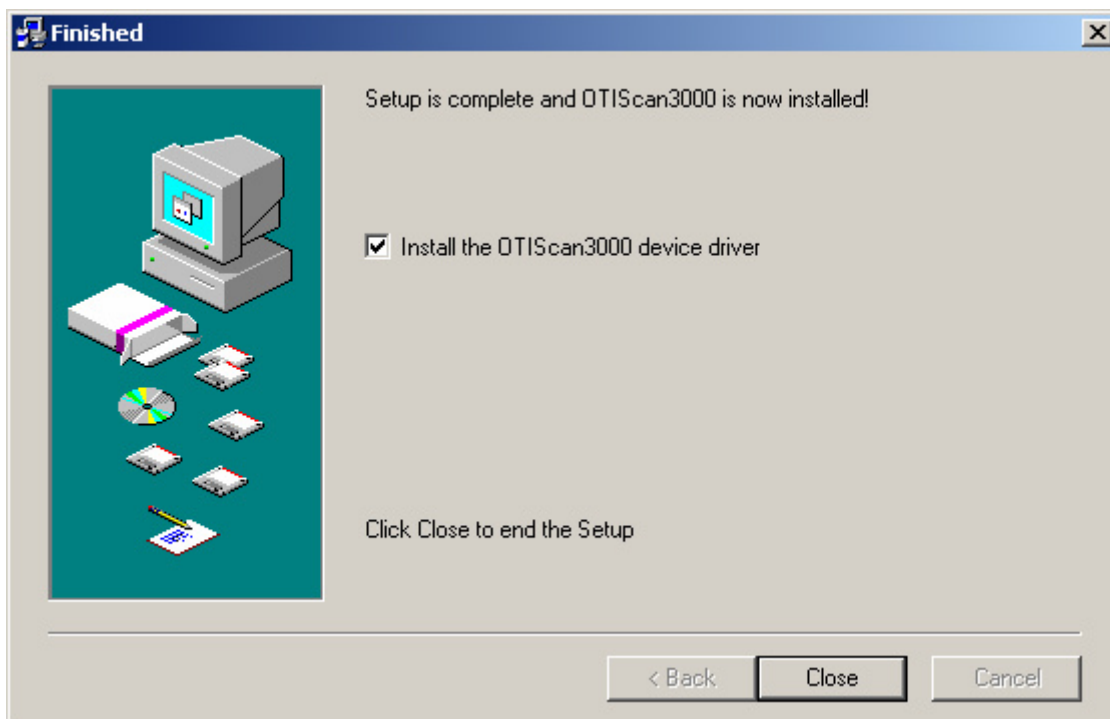


OTIScan 3000 Driver Installation


6. If you have chosen to install the OTIScan 3000 device driver, the Welcome screen of the VCP_1394 Driver installer will be shown. Make sure your OTIScan 3000 console is not connected to your computer, then click the Next button.

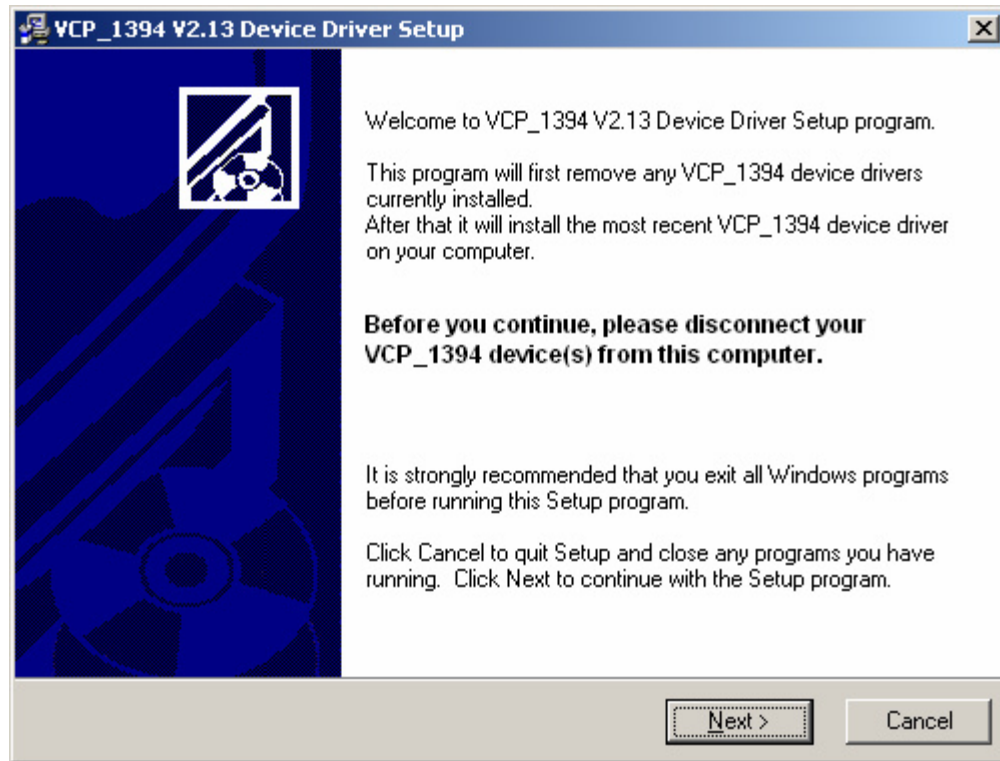

7. Click on the Next button on the Preparing to Install screen.

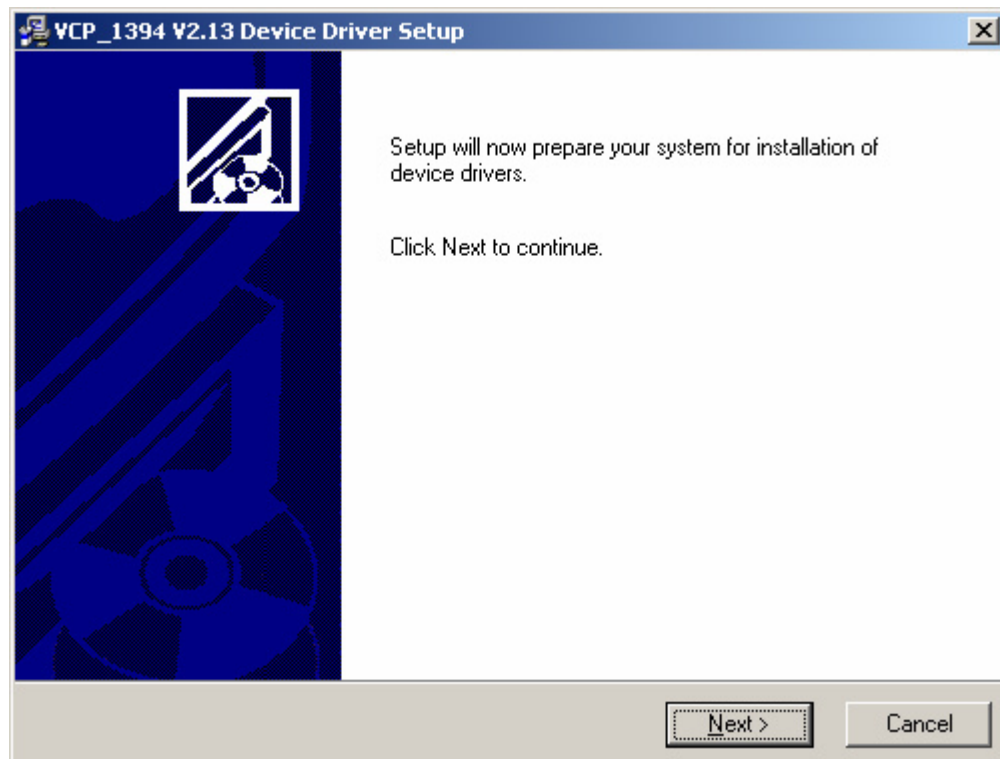


7. Click on the Finish button to exit the installer.

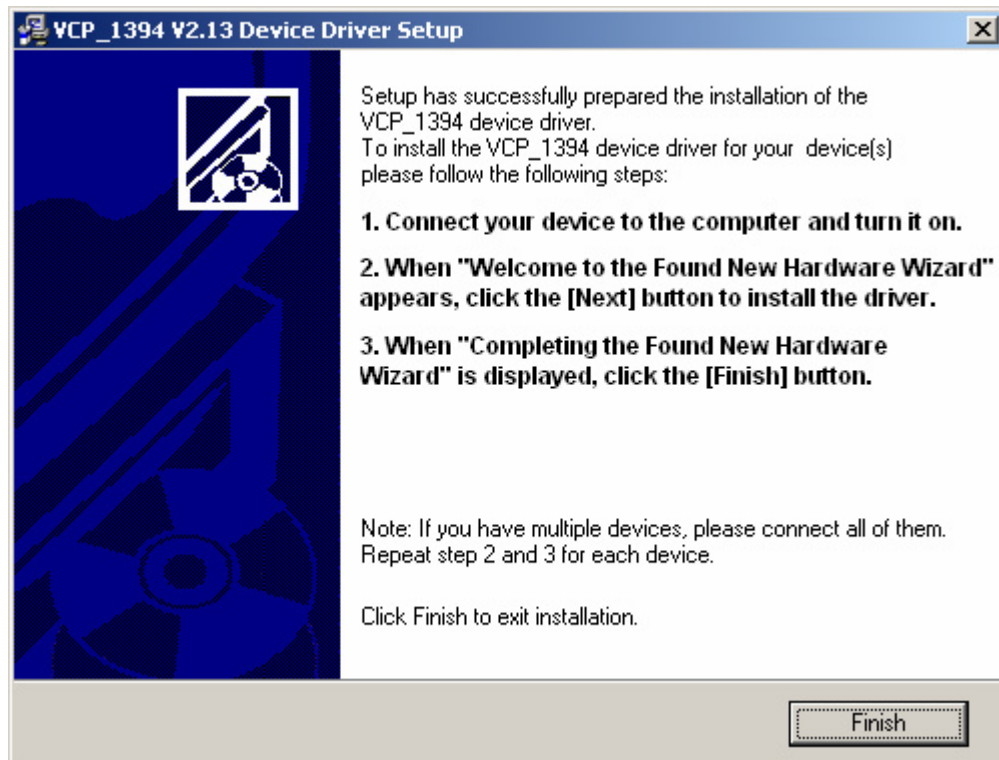


8. Connect your OTIScan 3000 console to your computer. On Windows Vista systems, the OTIScan 3000 driver installation is complete. If you are running Windows XP, the Found New Hardware Wizard will appear. Select "No, Not this time" and click the Next button.




9. Select "Install from a list or specific location" and click Next.

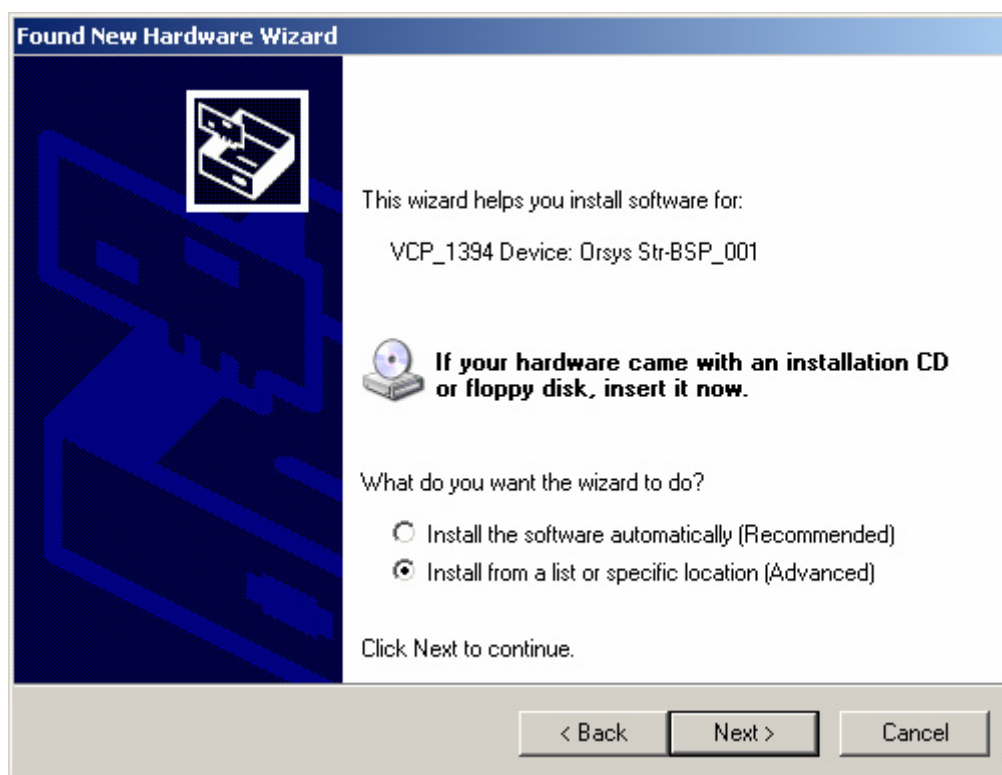

10. Select "Don't Search, I will choose the driver to install." and click Next.

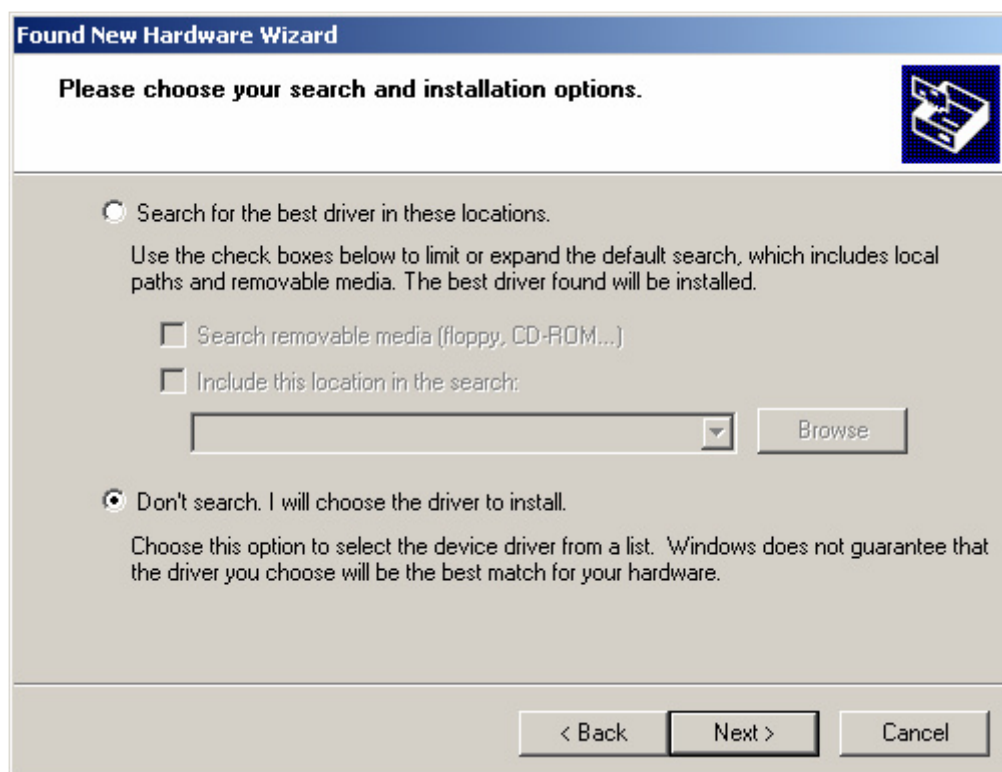

11. Select "VCP_1394 Device: Orsys Str-BSP_001" from the list and click Next.

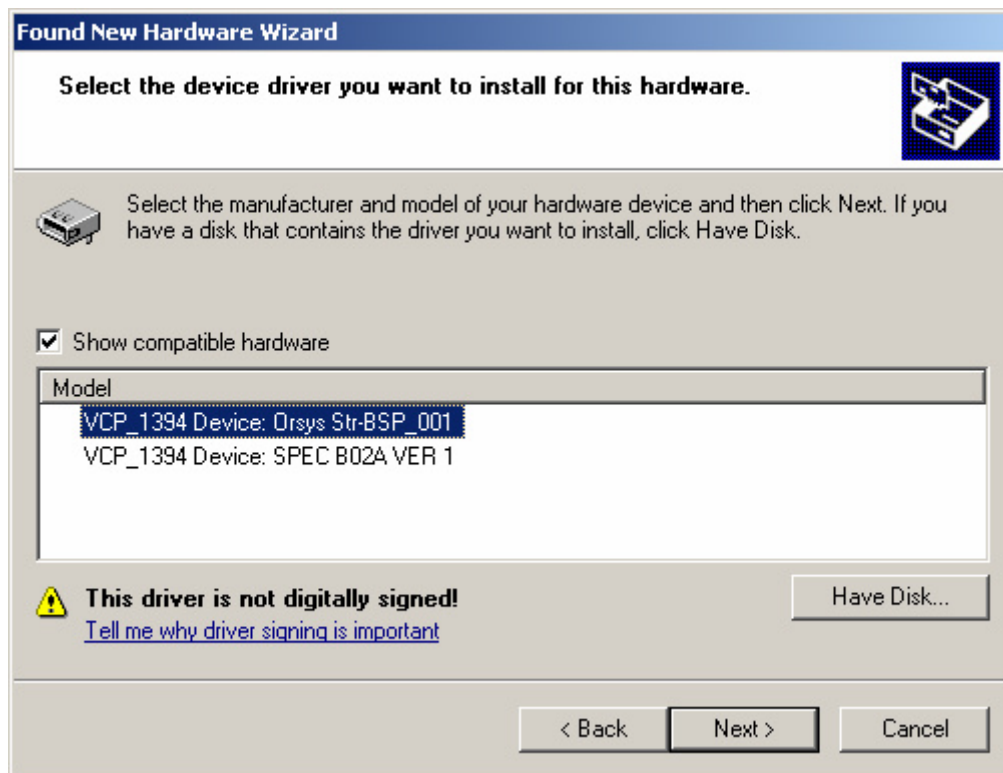


12. After the driver has been installed, click on the Finish button to close the Found New Hardware Wizard.



AScan Driver Installation

13. If you have received an AScan console, install the AScan drivers by double clicking on the RegisterAScanDrivers.bat file in the AScanDrivers directory on the CD.

14. You will be asked to confirm that your AScan console is not connected to your computer, then you will be warned that the driver has not undergone Windows Logo testing. Click "Continue Anyway." These dialogs will be shown twice.

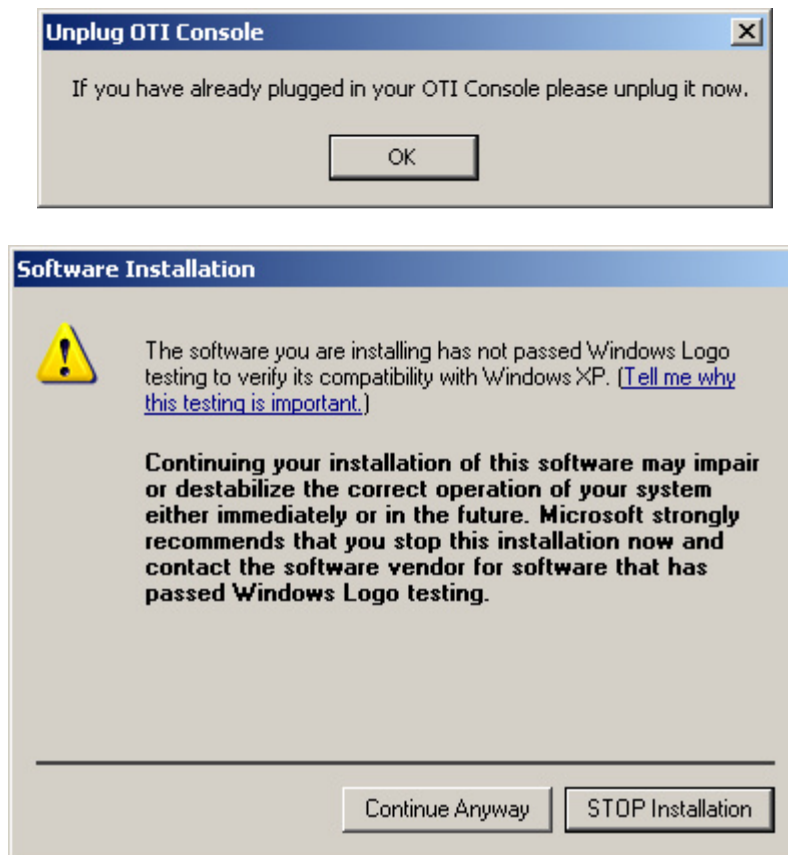

15. Connect your AScan console and turn it on. If you are running Windows Vista, AScan driver installation is complete. If you are running Windows XP, continue with the next step.


16. If you are running Windows XP, the Found New Hardware Wizard will appear. Select "No, Not this time" and click the Next button.




17. Select "Install from a list or specific location" and click Next.



18. Select "Don't Search, I will choose the driver to install." and click Next.

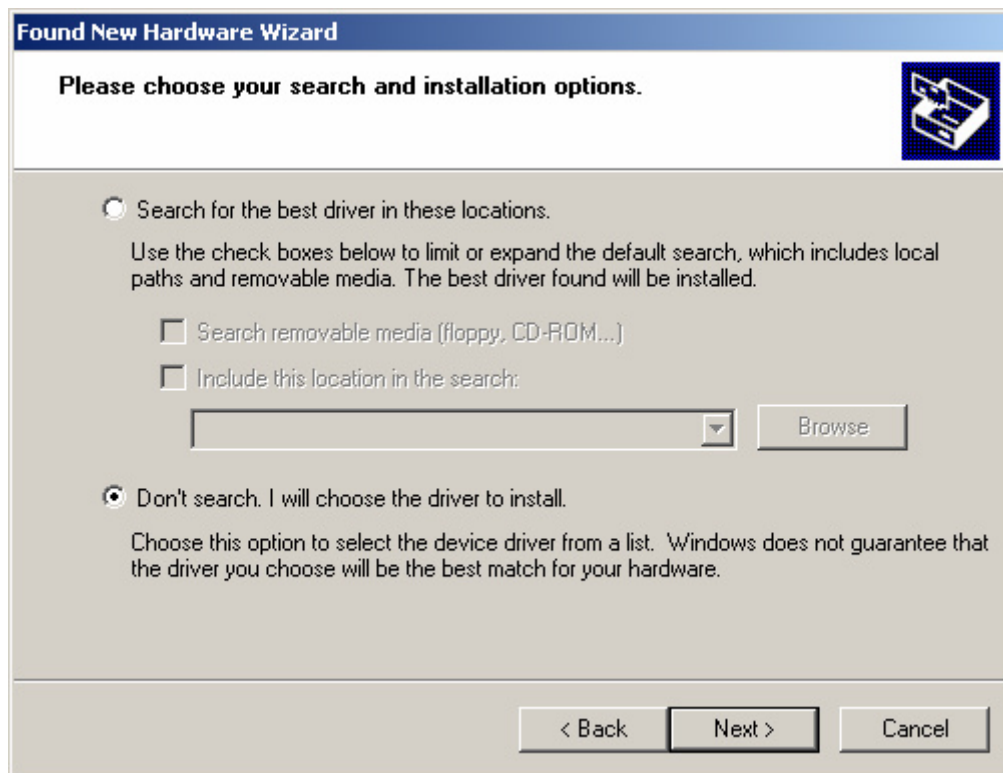

19. Select "USB High Speed Serial Converter" from the list and click Next.

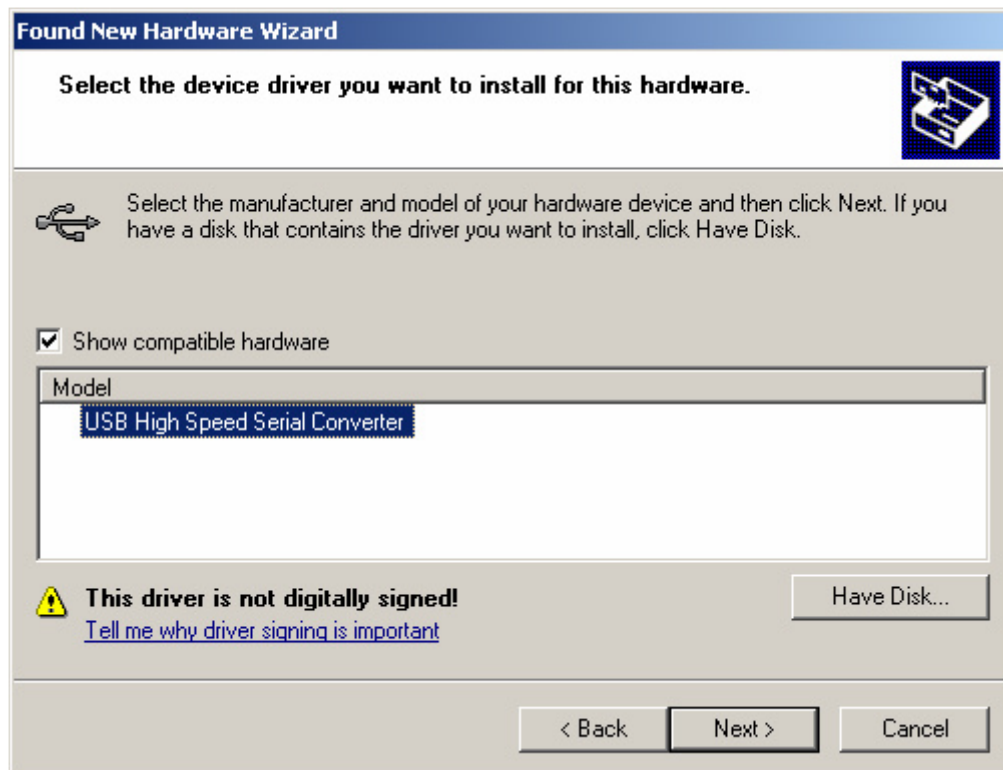


20. You will be warned that the driver has not undergone Windows Logo testing. Click "Continue Anyway."



21. After the driver has been installed, click on the Finish button to close the Found New Hardware Wizard.



22. The Found New Hard Wizard will appear again. Repeat steps 16 – 21 to install the “USB Serial Port” driver.

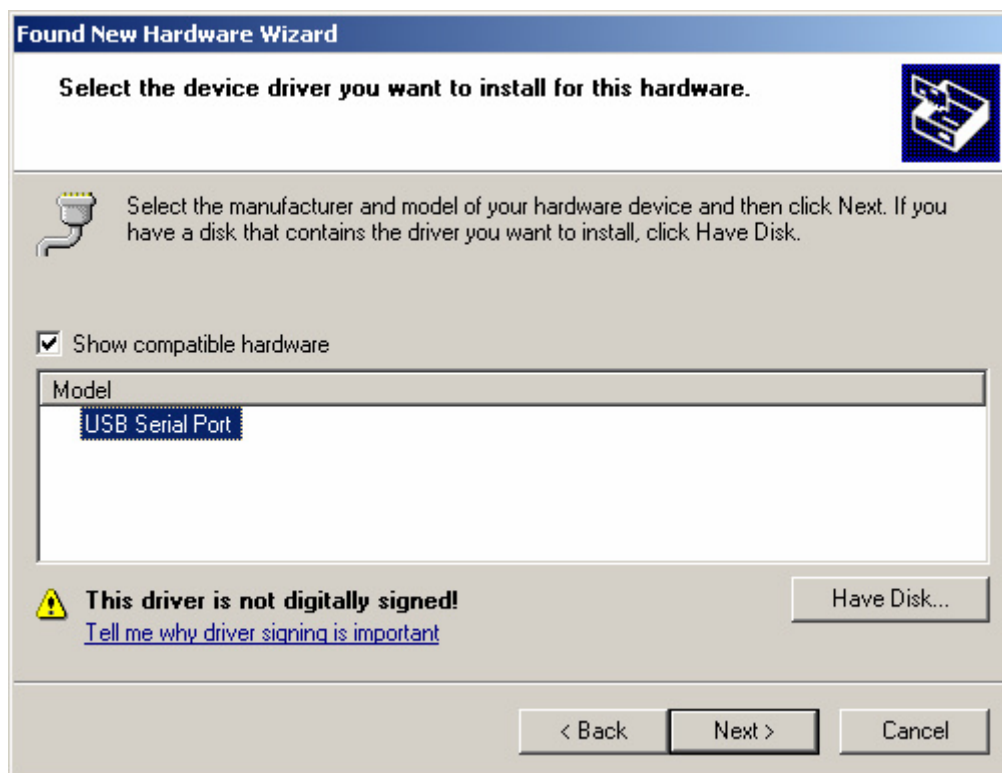

Copy msvc71.dll from CD

23. Run OTIScan 3000 by clicking on the shortcut in the Quick Launch bar or in the Start Menu under Programs | OTI | OTIScan 3000. If you receive an error message stating that msvc71.dll cannot be found, copy the file from the OTIScan 3000 CD to the C:\Program Files\OTI\OTIScan3000 folder.

OTIScan 3000 Uninstallation

1. To uninstall the OTIScan 3000 driver, run vcp_drv_uninstall.exe from the C:\Program Files\OTI\OTIScan3000\driver folder.
2. To uninstall the AScan drivers, run *Add or Remove Programs* from the Control Panel, select FTDI USB Serial Converter Drivers from the list, and click on the Change/Remove button.
3. To uninstall the OTIScan 3000 program, run *Add or Remove Programs* from the Control Panel, select OTIScan 3000 from the list, and click on the Change/Remove button.
4. After uninstalling OTIScan 3000, the following directories will remain on your hard drive. You may wish to delete them. Doing so will remove all patient data from your computer.

- C:\Library
- C:\Log
- C:\oti
- C:\Patient Log
- C:\Program Files\OTI\OTIScan3000

On Windows Vista you may also wish to delete the Compatability Files folder at:

C:\Users\<username>\AppData\Local\VirtualStore\Program Files\OTI\OTIScan3000